
Find out if Optune Gio® is right for you

How can Optune Gio® fit into my treatment plan if my doctor decides to prescribe Optune Gio®?

If you have newly diagnosed GBM, Optune Gio® is used together with chemotherapy* after surgery and radiation with chemotherapy.

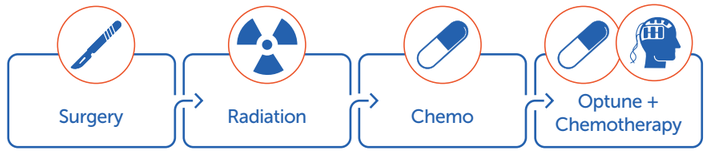

If you have recurrent GBM, Optune Gio® is used alone once treatment options like surgery and radiation have been exhausted.

Optune Gio® is not for everyone. For more information about who should not take Optune Gio®, please see the Optune Gio® User Manual.

How can Optune Gio® fit into my daily life once it has been prescribed?

Continuous treatment with Optune Gio® can be received almost anywhere. If your doctor says it's okay, you can use Optune Gio® during your daily routine at work, while travelling, and even at social functions.

*Temozolomide, also known as Temodal.®

Support throughout your journey

From your first question about Optune Gio® and throughout your treatment, there are many resources to help, including:

MyNovocure* support services provided by Novocure

Information on starting and staying on Optune Gio®

Visit **Optune.ca** for more information


*Please note, MyNovocure is formerly known as nCompass. The same support that healthcare professionals and patients know and trust is not changing, just evolving.

Patient images reflect the health status of the patients at the time each photo was taken.

novocure®

©2025 Novocure GmbH. All rights reserved.
Novocure, nCompass, MyNovocure, and Optune Gio® are registered trademarks of Novocure GmbH. All other trademarks, registered or unregistered, are the property of their respective owners.
CA-OPG-00001 May 2025 v1.0

OPTUNE GIO®

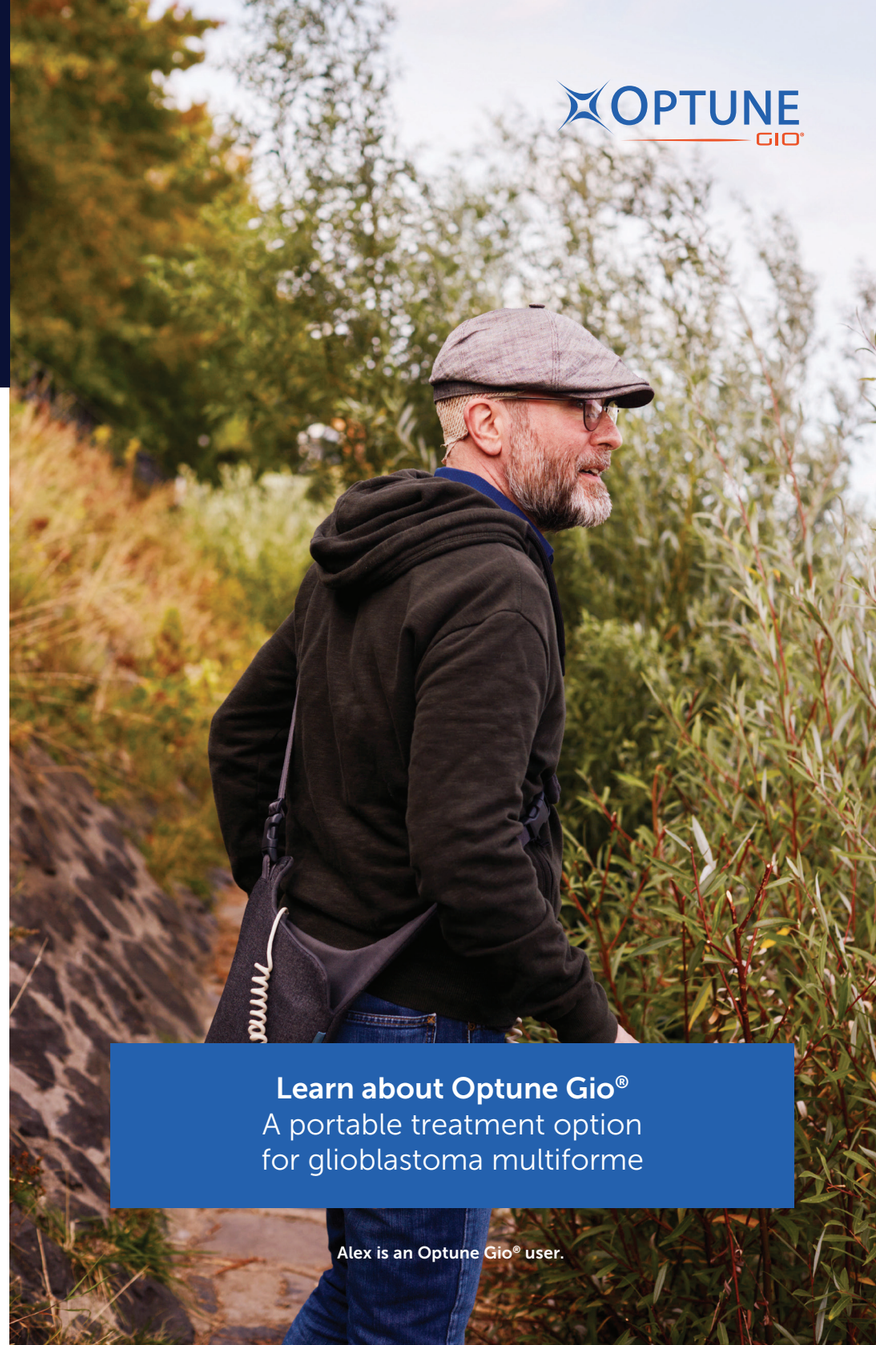

Learn about Optune Gio®
A portable treatment option for glioblastoma multiforme

Alex is an Optune Gio® user.

Learn the basics about Optune Gio®


What is Optune Gio®?

Optune Gio® is a portable medical device approved for the treatment of glioblastoma multiforme in adults aged 18 years of age or older.

How does Optune Gio® work?

Optune Gio® works by creating Tumour Treating Fields (TTFields), which are electric fields that are intended to prevent cancer cells from dividing.

TTFields can slow or stop GBM tumour cells from dividing and may cause some cells to die. Using 4 adhesive patches called transducer arrays, TTFields therapy is delivered to the brain where the cancer is located.

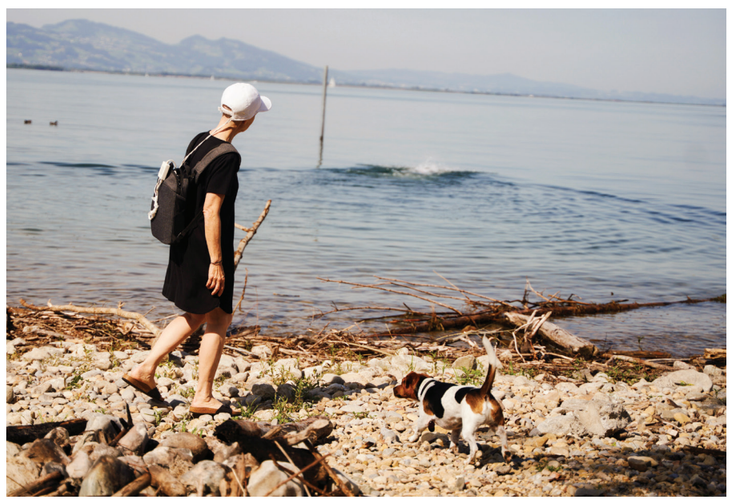

Andrea is an Optune Gio® user.

What is Optune Gio® approved to treat?

Optune Gio® is a portable medical device indicated to treat a type of brain cancer called glioblastoma multiforme in adult patients 18 years of age or older.

Newly diagnosed GBM

If you have newly diagnosed GBM, Optune Gio® is used together with and after standard of care maintenance chemotherapy if:

- Your cancer is confirmed by your healthcare professional AND
- You have had surgery to remove as much of the tumour as possible and completed radiation therapy

Recurrent GBM

If your tumour has come back, Optune Gio® can be used alone as an alternative to standard medical therapy if:

- You have tried surgery and radiation and they did not work or are no longer working AND
- You have tried chemotherapy and your GBM has been confirmed by your healthcare professional

Who should not use Optune Gio®?

Optune Gio® for GBM is not for everyone.
Talk to your doctor if you have:

- An implanted medical device (programmable shunt), skull defect (such as missing bone with no replacement), or bullet fragments. Optune Gio® has not been tested in people with implanted electronic devices, which may cause the devices not to work properly, and Optune Gio® has not been tested in people with skull defects or bullet fragments, which may cause Optune Gio® not to work properly
- A known sensitivity to conductive hydrogels (the gel on the arrays placed on the scalp like the ones used on ECGs). When Optune Gio® comes into contact with the skin, it may cause more redness and itching or may rarely cause a life-threatening allergic reaction

Do not use Optune Gio® for GBM if you are pregnant or are planning to become pregnant. It is not known if Optune Gio® is safe or effective during pregnancy.

What should I know before using Optune Gio®?

Optune Gio® should only be used after receiving training from qualified personnel, such as your doctor, a nurse, or other medical staff who have completed a training course given by Novocure®, the maker of Optune Gio®.

- Do not use any parts that did not come with the Optune Gio® Treatment Kit sent to you by Novocure or given to you by your doctor
- Do not get the device or transducer arrays wet
- Please be aware that Optune Gio® has a cord that may cause tripping when connected to an electric socket
- If you have an underlying serious skin condition on the scalp, discuss with your doctor whether this may prevent or temporarily interfere with Optune Gio® treatment

What are the possible side effects of Optune Gio®?

Common side effects of Optune Gio® when used together with chemotherapy for GBM were scalp irritation (redness and itchiness), low blood platelet count, nausea, vomiting, tiredness, and seizure.

The most common side effects when using Optune Gio® alone for GBM were scalp irritation (redness and itchiness) and headache. Other side effects were malaise, muscle twitching, fall, and skin ulcers.

Talk to your doctor if you have any of these side effects or questions.

Please see the Optune Gio® User Manual for complete information regarding the device's indications, contraindications, warnings, and precautions.

OPTUNE GIO®